

PRODUCT DESCRIPTION:

Each aliquot contains 10 µg of DNA extracted from a pure culture of *Haemophilus ducreyi*. The identification of this organism was confirmed by 16S sequencing. The purity of the culture was monitored by additional culturing and Gram staining to detect any contaminating bacteria. The DNA was extracted from the cells following a custom protocol utilizing Qiagen® Genomic DNA Handbook using Qiagen® Genomic DNA Buffers with a 500/G genomic tip. DNA concentration and A_{260/280} ratios are determined using a NanoDrop®. The extracted DNA also tested positive on an in-house real time PCR assay.

INTENDED USE:

Purified Genomic DNA is designed for use as an amplification and/or detection control for nucleic acid testing of *Haemophilus ducreyi*. It can also be used to determine a limit of detection (LOD), in diagnostic assay development, cross-reactivity studies or genomic sequencing. When used as a control for nucleic acid tests, the same protocols as those used to amplify extracted clinical specimens should be employed.

PRECAUTIONS:

- Use Universal Precautions when handling Genomic DNA.
- The material may be re-frozen after thawing. Repetitive freezing and thawing is not recommended (aliquot material if necessary).
- To avoid cross-contamination, use separate pipette tips for all reagents.

RECOMMENDED STORAGE:



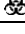
This control is supplied in TE Buffer and should be frozen at -20°C or below.

DO NOT USE IN HUMANS. FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

These products are intended for research, product development or manufacturing use only. These products are NOT intended for use in the manufacture or processing of injectable products subject to licensure under section 351 of the Public Health Service Act or for any other product intended for administration to humans.

© Registered trademarks are the property of their respective owners.

This product was manufactured in a facility which has a Quality Management System that is ISO 13485 certified.

REF	Catalog Number		Temperature Limitation
LOT	Lot Number		Expiration Date
RUO	For Research Use Only		Biological Risk

ZeptoMetrix Corporation
 878 Main Street
 Buffalo, NY 14202
 800-274-5487
 www.zeptometrix.com