



NATtrol™ *Chlamydia trachomatis* Positive Control Pack
Strain: LGV-II-434 Part Number: NATCT(434)-6MC-IVD

INTENDED USE:

NATtrol™ *Chlamydia trachomatis* (CT) Positive Control is an unassayed *in vitro* diagnostic external run control intended to be used with qualitative molecular assays for detection of nucleic acids from this organism. The control is intended to be used as an aid to diagnosis in that it is used to verify performance of the assays used to detect a physiological or pathological state. The routine and repetitive use of external run controls enables laboratories to monitor daily test variation, lot-to-lot test kit performance, individual operator variation, and can provide assistance in identifying increases in random or systemic errors. NATtrol™ *Chlamydia trachomatis* (CT) Positive Control contains intact organisms and should be run in a manner identical to that used for clinical specimens. This qualitative control is not automated and does not have an assigned value and it is the responsibility of the end user to establish their own target specifications for the control using their laboratory's molecular procedures.

PRODUCT SUMMARY AND EXPLANATION:

NATtrol™ CT Positive Controls are formulated with purified, intact organisms that have been chemically modified to render them non-infectious and refrigerator stable*.

Each NATtrol™ CT Positive Control Pack contains 6 x 1.25 mL vials of NATtrol™ *Chlamydia trachomatis* formulated in a Purified Protein Matrix that is fully commutable with true clinical specimens.

*Pat.: <http://www.zeptometrix.com/patent-information/>

PRINCIPLE:

The CT Positive Controls contain *Chlamydia trachomatis* cells inactivated by ZeptoMetrix's patented NATtrol™ process formulated in a proprietary Purified Protein Matrix that mimics the composition of a true clinical specimen. These are full process controls designed to monitor the effectiveness of extraction, amplification, and detection in nucleic acid testing procedures. These controls are suitable for use in in-house molecular assays and commercially available platforms.

WARNINGS AND PRECAUTIONS:

Although the NATtrol™ CT Positive Controls contain inactivated microorganisms, handling and disposal should be conducted as if potentially infectious.

This control contains material of human and animal origin and the user should observe Universal Precautions when handling and disposing of this product. Disposal must follow local regulations, if more stringent than regulations enforced by the CDC or the FDA.

Do not pipette by mouth.

To avoid cross-contamination, use separate transfer pipettes or tips for all materials.

Do not use beyond the expiration date shown on the label.

Usage beyond the expiration date shown on the label or after storage outside the recommended temperature may be detrimental to product performance or stability and lead to invalid or erroneous results.

If product is received damaged or leaking, contact ZeptoMetrix for instructions.

Changes in physical appearance of the product such as excessive turbidity, presence of precipitates, or discoloration may indicate degradation or contamination of the product. Discard the vial.

Failure to follow the assay or kit manufacturer's instructions explicitly for testing and analysis of results may lead to invalid or erroneous results.

If the expected result is not obtained, contact ZeptoMetrix for instructions.

NOT FOR USE IN HUMANS:

These products are NOT intended for use in the manufacture or processing of injectable products subject to licensure under the USA Food and Drug Administration Section 351 of the Public Health Service Act or for any other product intended for administration to humans.

RECOMMENDED STORAGE:

NATtrol™ CT Positive Controls should be stored at 2-8°C upon arrival.

When stored as directed, controls are suitable for use for up to 56 days (8 weeks) once opened.

INSTRUCTIONS FOR USE:

Vortex NATtrol™ CT Positive Control vials for 10 seconds to mix.

Follow the manufacturer instructions for use as a clinical sample.

Any serious incident that has occurred in relation to the device shall be reported to ZeptoMetrix and the competent authority of the Member State in which the user and/or the patient is established.

LIMITATIONS:

NATtrol™ CT Positive Controls are USA FDA Class 1 exempt, unassayed, *in vitro* diagnostic external run controls and are for professional use only.

NATtrol™ CT Positive Controls are not intended for use as a substitute for the internal controls provided by *in vitro* diagnostic kit manufacturers.

NATtrol™ CT Positive Control is not intended for use as a primary reference standard or material for any assay or testing procedure.

Quality control materials should be used in accordance with local, state, federal and accreditation requirements.

EXPECTED RESULTS:

One lot of NATtrol™ CT Positive Control was tested in two separate labs on a total of three different assays. Qualitative results are shown in Table 1.

As stated in the intended use, this product does not have an assigned value. Each laboratory must evaluate each lot of controls and establish acceptance criteria with their own specific molecular assay procedure and according to their own established quality assurance requirements and guidelines.

Product homogeneity has been demonstrated by validation studies and quality control testing.

Table 1:

Assay	Site	CT Results n=2	
		Urine	Swab
Cepheid Xpert® CT/NG	1	Positive	Positive
	2	N/T*	N/T*
BD ProbeTec™ ET CT/NG	1	N/T*	N/T*
	2	Positive	Positive
Hologic APTIMA Combo 2® CT/NG	1	N/T*	N/T*
	2	Positive	Positive

*Not Tested



NATtrol™ Chlamydia trachomatis Positive Control Pack
Strain: LGV-II-434 Part Number: NATCT(434)-6MC-IVD

ETIOLOGIC STATUS/BIOHAZARD TESTING:

NATtrol™ inactivation was completed on the stocks used to formulate each control and further verified by the absence of bacterial growth in a validated growth protocol.

The Purified Protein Matrix used in the manufacture of this product is treated with 0.09% sodium azide. It was manufactured from materials that have been tested and found non-reactive at the donor level for HIV-1/HIV-2 Antibody, HBsAg and HCV Antibody by FDA licensed donor screening test methods. All materials are also tested for HIV-1 and HCV by FDA approved Nucleic Acid Test (NAT) methods. Heat inactivated bovine based source materials used in the manufacture of this product meet applicable USDA requirements for abattoir sourced animals, traceability and country of origin. The materials were collected at USDA licensed establishments or legally imported from countries recognized by the USDA as negligible or controlled for risk for Bovine Spongiform Encephalopathy (BSE) and other exotic disease agents. Donor animals were inspected ante and post mortem at the abattoir as required by the USDA.

PRODUCT WARRANTY:

ZeptoMetrix LLC's limited product warranty and other terms and conditions related to the purchase and use of ZeptoMetrix products are set forth in ZeptoMetrix's Terms and Conditions of sale found on ZeptoMetrix's website at [Sales Terms and Conditions](#). If you have any questions, please contact ZeptoMetrix Customer Service at zepto.customerservice@antylia.com.

DISCLAIMER AND LIMITATION OF LIABILITY:


ZeptoMetrix LLC disclaims all warranties with respect to this document and the information contained herein, expressed or implied, including but not limited to those of merchantability, fitness for a particular purpose, or non-infringement. To the extent allowed by law, in no event shall ZeptoMetrix LLC be liable, whether in contract, tort, warranty, or consequential damages or lost profits in connection with or arising from this document and the information contained herein, including but not limited to the use thereof even if ZeptoMetrix is advised of the possibility of such damages.

LEGEND OF LABELING SYMBOLS:

	Manufacturer		Temperature Limitation
	In vitro Diagnostic Use		Use-By Date
	European Mark of Conformity		Biological Risk
	Part / Catalogue Number		Authorized Representative
	Batch Code		Consult Instructions for Use
	Contains biological material of animal origin		Contains biological material of human origin
	Unique Device Identifier		

Manufacturer:
 ZeptoMetrix LLC
 25 Kenwood Circle
 Franklin, MA 02038, USA

EC Representative:

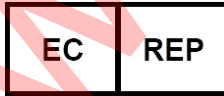

EMERGO EUROPE
 Westervoortsedijk 60
 6827 AT, Arnhem
 The Netherlands

UK Responsible Person:
Emergo Consulting (UK) Limited
 c/o Cr360 – UL International,
 Compass House, Vision Park Histon,
 Cambridge CB24 9BZ, United Kingdom

REVISION HISTORY.

Revision Level	Description of Revisions
11	Added additional statement to warning concerning usage beyond the expiration date shown on the label or after storage outside the recommended temperature, changes in physical appearance of the product, failure to follow the assay or kit manufacturer's instructions explicitly, and if the expected result is not obtained action. Added "NATtrol™ CT Positive Control is not intended for use as a primary reference standard or material for any assay or testing procedure." to Limitations. Added "Any serious incident that has occurred in relation to the device shall be reported to ZeptoMetrix and the competent authority of the Member State in which the user and/or the patient is established." to Instructions for Use. Added verbiage to Intended Use to indicate it is a qualitative control that is not automated. The fundamental Intended Use has not changed. Added Revision History section. Changes in response to BSI Technical File review for IVDR compliance. Added biological and UDI symbols. Updated EC Representative Emergo Europe address to new location. Added UK Responsible Person contact details.